


Cat. No	Description	Articles
---------	-------------	----------

<Continued from... 50 Lit. Digital Rotary Evaporator, "WEV-1050" >

**DAIHAN® Digital Rotary Evaporator Full System**

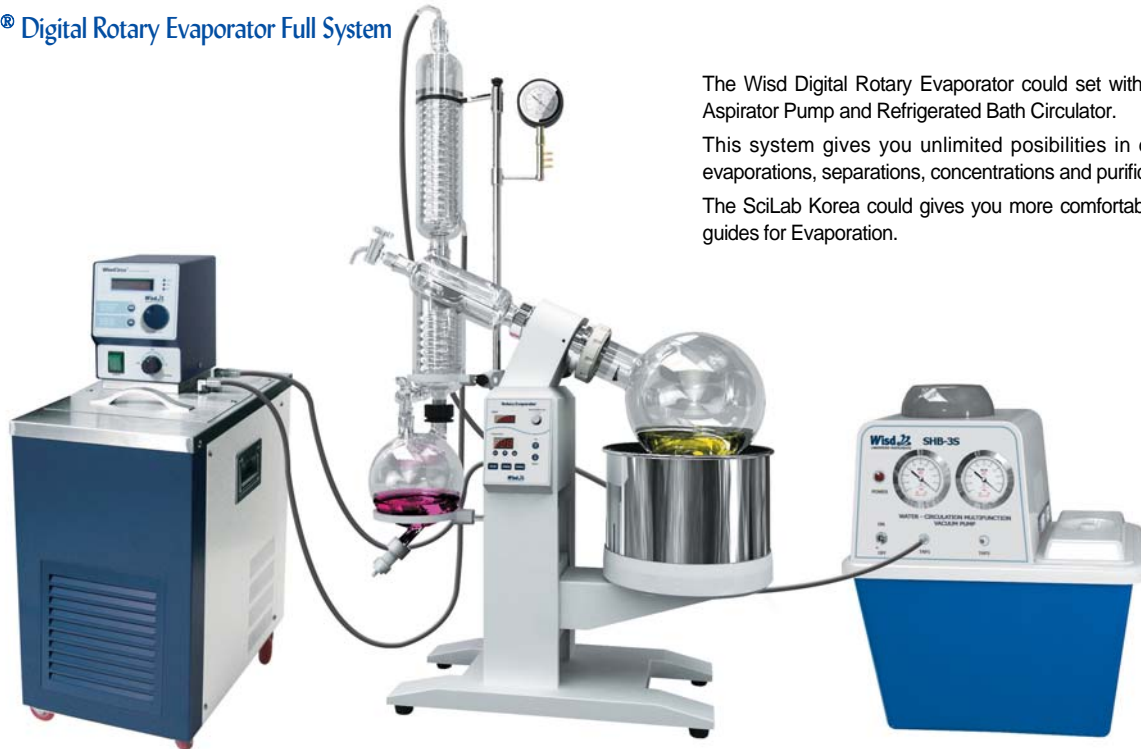

The Wisd Digital Rotary Evaporator could set with Circulation Aspirator Pump and Refrigerated Bath Circulator. This system gives you unlimited possibilities in distillations, evaporations, separations, concentrations and purifications. The SciLab Korea could gives you more comfortable selection guides for Evaporation.

Evaporator Full System : Digital Rotary Evaporator + Circulation Aspirator Pump + Refrigerated Bath Circulator  
 Each Items Details See Page P. XXX, P. XXXX

**DAIHAN® Electric Circulation Aspirator Pump, "WEV-003S" 15 Lit. PP Water Tank, Safety Check Valve NEW**  
 80 Lit./min, 0.098 MPa



**Wisd**  
 2-Year Warranty

- DAIHAN® Designed China-made
- Cyclically use Water, Save the Precious Water Resources.
- Save more than 35% of Electricity than other type of Vacuum Pump.
- Double-tap, Double-display, can be used alone or in Parallel.
- Resist Corrosion by Acid, Alkali and Solvents.

- \* The Fluid Muffler Purposemade can Reduce the Gas in Water and make the Vacuum Degree Higher and more Stable Reduce the Friction between Gas and Fluid, Reduce the Noise.
- \* Iuorine Rubber Sealing, the inner of which can't be Intruded by Corrosive Gas.
- \* Need to Replace the Water in the Flume Regularly to ensure the Purity of Water Quality, the Vacuum Degree and to Avoid Dirt Stains.
- \* Can be used to Extract Corrosive Gas, need to Shorten the Period of Water Changing.

**Specification**

Spec.	Model	"WEV-003S"
Capacity(Water Tank)		15 Lit.
Max. Vacuum Degree		0.098 MPa
Flow		80 Lit./min
Lift		10 m
Single Tap Air Sucking Amount		10 Lit./min, 2A
Safety Function		Check Valve
Material of Water Tank		PP
Dimension & Net Weight		385 × 280 × h420 mm / 11 kg
Power Consumption		180 W
Power Supply		220V, 50Hz

**Electric Circulation Aspirator Pump "WEV-003S"**  
 15 Lit., 0.098MPa, 80Lit./min

**Ordering Information**

DH.WEV0003S Circulation Aspirator Pump, 15Lit., "WEV-003S", 220V